





### O-6 WRIST SUPPORT

8 cm width special woven elastic band, one end is sewed in a triangular shape with a velcro providing required compress. 5 Sizes are available.

**Indication:** Used for wrist and traumas, sprain and rheumatic pains, sports activities to protect the wrist from injuries.



### O-11 WRIST ORTHOSIS

Made from an elastic fabric, it can be adjusted from 4 different points by velcros, has specially shaped aluminium support at the bottom part of orthosis affixed into pocket to keep the hand in resting position. Supports can be shaped easily according to the expert. 5 sizes are available. One does not fit both hands.

**Indication:** Used for wrist traumas, sprains, injuries, and keeping the hand at the resting position after operations.



### O-11A WRIST ORTHOSIS

Orthosis with a thumb loop is designed in accordance with wrist's anatomic shape in order to provide a precise fitting. It doesn't fit both hands. 5 sizes are available.

**Indication:** Used on wrist traumas, sprains, injuries, and keeping the hand in the resting position after operations.



### O-11B WRIST THUMB ORTHOSIS

In addition to the orthosis which keeps the hand in resting position, it has a thumb support immobilizing the thumb. The stabilization is provided by shaping the aluminium support beside the thumb. It doesn't fit both hands. 5 sizes are available.

**Indication:** Used when the thumb should be relaxed.

### O-11C WRIST ORTHOSIS WITH BANDAGE

An additional support is provided to the orthosis, which is made from specially woven elastic skin color fabric, by an elastic band. Different models for right and left hand and 5 different sizes are available.

**Indication:** Used on wrist traumas, sprains, injuries, and keeping the hand in the resting position after operations.



### O-11R RIGID WRIST ORTHOSIS

Made from a special fabric, the orthosis is designed in accordance with the wrist anatomy. Closure is provided by velcros from the three points of the orthosis. The stabilization is supported by the velcros which grasp around the plastic buckles and make puller. The extra support is provided by the 2 plastic baleens putted on the orthosis. An aluminium support which is shaped by the molds that are put into the pocket. It immobilizes the wrist in neutral resting position. 4 sizes are available in dark blue and skin colors. Also right and left models are available.

**Indication:** Used for wrist injuries, muscles and ligament pains, Carpel Tunnel illnesses after operations.



### K-11 AIRY WRIST ORTHOSIS

The airy wrist orthosis, suitable for the wrist anatomy provides closure from 2 different points by velcros. It is made from special perforated padded fabric. Additionally the velcro band entirely wrapping around the wrist provides an extra support. Velcros go through the plastic buckles and turns back from them, pulls the orthosis and provides a precise guidance to the wrist. Standard aluminium splints are inserted into the pocket, keeps the wrist in neutral resting position and immobilizes it. 3 sizes as right and left models are available in black color.

**Indication:** Used for wrist sprains and pains, ruptures of muscle and ligaments, carpal tunnel syndrome and after operations.

### K-11A AIRY WRIST THUMB ORTHOSIS

The thumb orthosis, suitable for the wrist anatomy, is designed from special perforated padded fabric as immobilizing the wrist and thumb, covers the thumb as well. The aluminium splint at the underpart keeps the wrist in resting position and another splint by one side keeps the fingers in resting position. 3 sizes as right and left models are available in black color.

**Indication:** Used for wrist and thumb muscle and ligament twistings, de querven tendinitis, carpal tunnel syndrome, weak muscles and spasticity in order to prevent wrist and thumb contractures.



### M-7 ELASTIC WRIST ORTHOSIS

Made from an elastic, air-permeable fabric prevents the body from perspiring, it can be adjusted from 4 different points by velcros, has specially shaped aluminium support at the bottom part of orthosis affixed into pocket to keep the hand in resting position. Supports can be shaped easily according to the expert. 5 sizes are available. One does not fit both hands.

**Indication:** Used for wrist traumas, sprains, injuries, and keeping the hand at the resting position after operations.



### M-8 ELASTIC WRIST SUPPORT

Made from an elastic, air-permeable 8 cm width special fabric prevents the body from perspiring. One end is sewed in a triangular shape with a welcro providing required compress. 5 Sizes are available.

**Indication:** Used for wrist and traumas, sprain and rheumatic pains, sports activities to protect the wrist from injuries.

### O-12 FOREARM ORTHOSIS

The forearm orthosis made from a 27cm-thick elastic cloth is about 10cm longer than wrist orthosis and spreads over the forearm part as well. One model with thumb loop and one without thumb loop and 5 sizes are available. It doesn't fit both right and left hand.

**Indication:** Used on forearm and wrist traumas and used when forearm should be relaxed.



### O-13 DYNAMIC APPARATUS

The apparatus made from 2mm thick hard aluminium has wrist joint. Hard elastic band engaged to the wrist joint gets the wrist exercised. 4 steely strings are affixed onto hand. Strings have shafts at the fingers' tip. The steely string which is affixed to the thumb comes from one side. Steely strings can easily be flexed and adapted to fingers. Inner and outer parts are lined with 3 mm plastizote. Right hand and left hand models are available for man size and woman size.



**Indication:** Used after surgical operations; on paralytic patients and physio-therapies to strengthen the finger muscles.



### O-14 STATIC APPARATUS

Aluminium and vitraten two different models are available. Inside is lined with 3mm plastizote. It can be contacted from three different points by velcros. Right hand and left hand models are available for man size and woman size.

**Indication:** Used on paralytic patients; for wrist illnesses; after plaster application; when the hand should be kept at the resting position.

### O-21 ALUMINIUM APPARATUS FOR FINGER

1mm thick 2cm width, edges of the apparatus are rounded aluminium orthosis. Inside is lined with plastizote.

**Indication:** Used to immobilize the finger after fractures and deformities.



### O-22 ALUMINIUM APPARATUS FOR FINGER "T"

1mm thick 2cm width, edges of the apparatus are rounded aluminium orthosis. Inside is lined with plastizote. "T" orthosis has an additional immobilizer and produced 33 cm length.

**Indication:** Used to immobilize the finger after fractures and deformities.



### O-23 ALUMINIUM APPARATUS FOR FOINGER WITH WRIST SUPPORT


It is designed with an additional brace with a special pocket where the apparatus can be affixed in order to provide an ease of immobilization

**Indication:** Used for finger fractures, ekstansör tendon ruptures, and finger deformities

### O-24 FINGER ORTHOSIS

It is produced in special plastic molds. 8 sizes are available. It has small scuttles and nail parts are open.

**Indication:** Used for finger fractures, ekstansor tendon ruptures, and finger deformities.

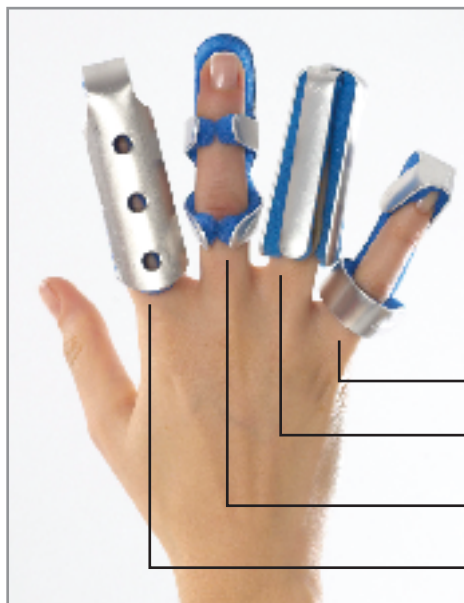

### O-41 APPARATUS FOR FINGER

Made of aluminium inside part of the pressed apparatus is lined with sponge. It can be immobilized without using any band easily. Each apparatus has 3 different sizes.

**Indication:** Used on finger fractures and deformities, sports injuries to immobilize the finger

Baseball Splinti

Quadruplet Splint

Frog Splint

Binary Splint



### O-50 DYNAMIC APPARATUS FOR FINGER

Made from stainless steel spring, the points which contact the finger is rounded by special plastics. The middle point of the spring is mounted onto the plastics with a screw and the other points are mounted onto the plastics with stainless rivets.

**Indication:** Used for tendon and nerve injuries.



### N-11 WRIST/THUMB SUPPORT

Orthosis made of specially woven elastic skin color fabric is designed in accordance with wrist/thumb anatomy. It has a real leather pocket in the palm, and a reshapable aluminium support is inserted. A heart shape pocket sewed on thumb part with a plastic support inside provides stability. Different models for right and left hand and 5 different sizes are available.

**Indication:** Used after bandage application and surgical operations; for Dekuarven and Carpel tunnels illnesses; and when the thumb and wrist should be relaxed.



### X-5 THUMB OTHOSIS

The orthosis made of elastic fabric lined with combed cotton is designed according to thumb's anatomic shape. A shapable soft aluminium support can be affixed into small pockets nearby thumb and it can be immobilized with a welcro on the wrist. Dark-blue, 5 sizes are available.

**Indication:** Used for thumb ligament injuries and ruptures, when thumb should not move.



### B-11 WOOLEN ELASTIC WRIST SUPPORT

It is made of elastic cloth and keeps the wrist warm and assist the treatment

**Indication:** Used for wrist and traumas, sprain and rheumatic pains, sports activities to protect the wrist from injuries