



K-1 WONROSSEN APPARATUS


The adjustment mechanization of the apparatus, is provided by velcros sewed on polyethilen belts in the legs and waist, and waist and femoral bands of the mechanization of the apparatus are metal. Polyethilen's edges are lined with soft vinyl. It can be applied on diaper and cleaned easily.

Indication: Used on dislocation of waist by doctor recommendation.

K-2 DENNIS BROUN APPARATUS

In the models with boot, an open finger boot, which is assembled to the apparatus, with puller and laces is used. The boot that finger part can be opened completely, provides availability to watch if child's foot fits precisely to the boot. The puller band is joint before boot is laced up and prevents foot to slip out the boot. The apparatus is made of aluminium as two separate rail pieces that one slides on another. Gear systems jointed to the ends of the pieces provide inside or outside rotation. Gear mechanization in the heel, jointed by screws and nails, makes it possible to assemble the shoe to the device with a required angle. Two models of the device, long and short one, are available. Boots are produced between 14 and 25 numbers.

Indication: Used on dislocated hip by doctor recommendation.

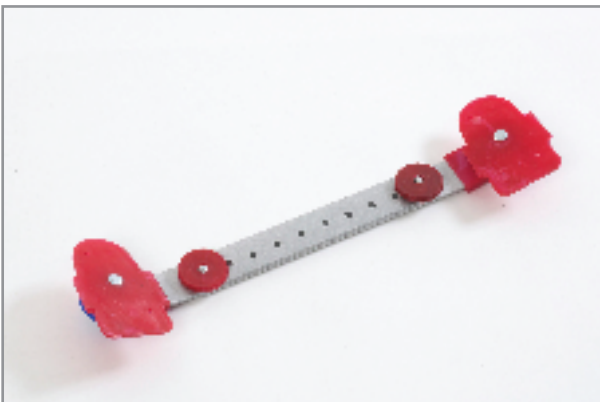

K-2A DENNIS BROUN APPARATUS WITHOUT BOOT

The apparatus is made of aluminium as two separate rail pieces that one slides on another. Gear systems jointed to the ends of the pieces provide inside or outside rotation. Gear mechanization in the heel, jointed by screws and nails, makes it possible to assemble the shoe to the device with a required angle. Two models of the device, long and short one, are available.

Indication: Used on dislocated hip by doctor recommendation.



K-3 PAVLIK BANDAGE

The parts on the shoulder and waist are made of a tightly woven cloth and inner part of that is lined with a soft material. The foot insertion section is sewed velcros so that they provide adjustment mechanization. The straps coming from shoulder and foot to the waist are easily adjusted by nettle buckles. 5 sizes are available.

Indication: Used on dislocated hip by doctor recommendation.

K-3A PAVLIK BANDAGE WITH VELCRO

The parts on the shoulder and waist are made of a padded fabric on which velcros can be stick and inner part of that is soft. The foot insertion section is sewed velcros so that they provide easy adjustment mechanization. 5 sizes are available.

Indication: Used on dislocated hip by doctor recommendation.



K-4 ABDUCTION CORSET

The edges of corset made from polyethilen, is lined with unvulcanized rubber and vinyl. It is immobilized by straps with velcro at the waist and shoulder, can be washed and wiped easily. 5 sizes are available.

Indication: Used for dislocated hip by doctor recommendation.



K-5 ABDUCTION PANT “FREJKA”

It is made by lining the inside of polyethilen with sponge and water-proof nylon. The corset is long enough in perianal region to abduct both two legs of infant. After perianal region, it is curved and opposite two ends are connected over the leg with velcro bindings. The edges of the corset which have air holes are strenghtened by sponge and lined with vinyl. 5 sizes are available.

Indication: Used for dislocated hip by doctor recommendation.



K-5A FREJKA APPARATUS OILCLOTH

The corset, which is produced from water-proof cloth, is put steel and aluminium balens inside and supported by sponge. It can be easily adjusted by buckles and straps. Femoral part has extra sponge. 5 sizes are available.

Indication: Used on dislocated hip by doctor recommendation.



K-7 ABDUCTION APPARATUS

Made from 15 mm metal that is nicked plated at back and front. There are enough holes to make the adjustment for the back and front parts of the apparatus. With the help of the tap holes, it is provided to find the correct place with screw to combine the parts of the apparatus. The inner side of the 3 mm plastic is covered with plastizote. All the edges of the plastic are sewed with vinylex by being put cords to the edges of the plastic on the calf part. The adjustment on the waist and calf is provided by the velcros. 5 sizes of the apparatus are available.

Indication: Used on dislocated hip by doctor recommendation. Usage is told by a specialist to the family.



X-28 ABDUCTION SHORT

The corset is made of 12 mm thick English made plastizote in skin color. Its all edges are rounded. And the neck support is preshaped by heat in the molds which are suitable for neck's anatomy. Skin can breath through holes on the plastizote. Diamond patterns on the plastizote provides a better stabilization for corset. At the waist, it is stabilized by 5-cm long velcro and at the femoral it is stabilized by two way 2-cm long velcro. 5 sizes ara available.

Indication: Used for dislocated hip by doctor recommendation.